

TECHNICAL DATA SHEET

# APC Anti-Mouse CD45 (30-F11)

Catalog Number: 20-0451

## PRODUCT INFORMATION

**Contents:** APC Anti-Mouse CD45 (30-F11)


**Isotype:** Rat IgG2b, kappa

**Concentration:** 0.2 mg/mL

**Clone:** 30-F11

**Reactivity:** Mouse

**Formulation:** 10 mM NaH<sub>2</sub>PO<sub>4</sub>, 150 mM NaCl, 0.09% NaN<sub>3</sub>, 0.1% gelatin, pH7.2

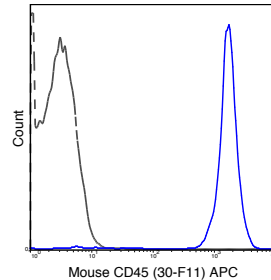

C57Bl/6 splenocytes were stained with 0.125 ug Anti-Mouse CD45 APC (20-0451) (solid line) or 0.125 ug Rat IgG2b APC isotype control (dashed line).

## DESCRIPTION

The 30-F11 antibody reacts with mouse CD45, which is one of the most abundant hematopoietic markers and is expressed on all leukocytes (the Leukocyte Common Antigen, LCA). CD45 is a protein tyrosine phosphatase existing in several isoforms, each being generated and expressed in cell-specific patterns. With its broad cell distribution, CD45 is critical for many leukocyte functions, regulating signal transduction and cell activation associated with the T cell receptor, B cell receptor, and IL-2 receptor. Other forms of CD45, with restricted cellular expression, include CD45R (B220), CD45RA, CD45RB, CD45RO and others. The 30-F11 antibody is widely used as a leukocyte marker for B cells, T cell subsets and NK cell subsets.

## PREPARATION & STORAGE

This monoclonal antibody was purified from tissue culture supernatant via affinity chromatography. The purified antibody was conjugated under optimal conditions, with unreacted dye removed from the preparation. It is recommended to store the product undiluted at 4°C, and protected from prolonged exposure to light. Do not freeze.

## APPLICATION NOTES

This antibody preparation has been quality-tested for flow cytometry using mouse spleen cells, or an appropriate cell type (where indicated). The amount of antibody required for optimal staining of a cell sample should be determined empirically in your system.

## REFERENCES

Panchal RG, Ulrich RL, Bradfute SB, Lane D, Ruthel G, Kenny TA, Iversen PL, Anderson AO, Gussio R, Raschke WC, and Bavari S. 2009. 284: 12874-12885. (Immunoprecipitation). Cherpes TL, Busch JL, Sheridan BS, Harvey SAK, and Hendricks RL. 2008. J. Immunol. 181: 969-975. (Complement-mediated cell depletion). Cheng G, Zhang H, Yang X, Tzima E, Ewalt KL, Schimmel P, and Faber JE. 2008. Am. J. Physiol. Regul. Integr. Comp. Physiol. 295: R1138-R1146. (Immunohistochemistry – paraffin embedded tissue). Nguyen JT, Evans DP, Galvan M, Pace KE, Leitenberg D, Bui TN, and Baum LG. 2001. J. Immunol. 167: 5697-5707. (Immunofluorescence microscopy, Immunoprecipitation). Czyzyk J, Leitenberg D, Taylor T, and Bottomly K. 2000. Mol. Cell. Biol. 20(23): 5740-5747. (Western Blot). Tsuboi S and Fukuda M. 1998. J. Biol. Chem. 273: 30680-30687. (Western Blot).