

TECHNICAL DATA SHEET

APC Anti-Human CD279 (PD-1) (J110)

Catalog Number: 20-2797

PRODUCT INFORMATION

Contents: APC Anti-Human CD279 (PD-1) (J110)

Isotype: Mouse IgG1

Concentration: 5µl (1µg)/test


Clone: J110

Reactivity: Human

Use By: 12 months from date of receipt

Storage Conditions: 2-8°C protected from light

Formulation: 10 mM NaH₂PO₄, 150 mM NaCl, 0.09% NaN₃, 0.1% gelatin, pH7.2

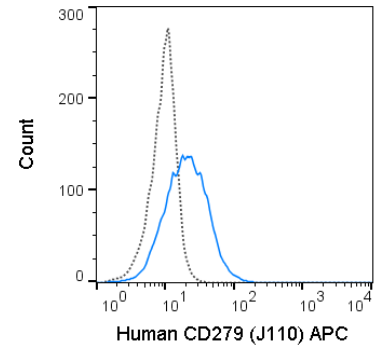

Human PBMC were stimulated for 3 days with PHA and stained with 5 µL (1 µg) APC Anti-Human CD279 (20-2797) (solid line) or 1 µg APC Mouse IgG1 isotype control (dashed line).

DESCRIPTION

The J110 antibody is specific for human CD279, also known as programmed death-1 (PD-1). PD-1 is a 55 kDa glycoprotein member of the Ig superfamily and is expressed on activated T cells, B cells, and myeloid cells. PD-L1 (B7-H1) and PD-L2 (B7-DC) are ligands for CD279 and belong to the B7 gene family. PD-1 is an immune checkpoint and interaction with its ligands results in inhibition of T cell proliferation and cytokine secretion.

PREPARATION & STORAGE

This monoclonal antibody was purified from tissue culture supernatant via affinity chromatography. The purified antibody was conjugated under optimal conditions, with unreacted dye removed from the preparation. It is recommended to store the product undiluted at 4°C, and protected from prolonged exposure to light. Do not freeze.

APPLICATION NOTES

This monoclonal antibody was purified from tissue culture supernatant via affinity chromatography. The purified antibody was conjugated under optimal conditions, with unreacted dye removed from the preparation. It is recommended to store the product undiluted at 4°C, and protected from prolonged exposure to light. Do not freeze.

REFERENCES

Freeman GJ, Long AJ, Iwai Y, et al. 2000. J Exp Med. 192:1027-1034. Latchman Y, Wood CR, Chernova T, et al. 2001. Nat Immunol. 2(3):261-268. Bennett F, Luxenberg D, Ling V, et al. 2003. J Immunol. 170(2):711-718.

Tonbo Biosciences tests all antibodies by flow cytometry. Citations are provided as a resource for additional applications that have not been validated by Tonbo Biosciences. Please choose the appropriate format for each application and consult Materials and Methods sections for additional details about the use of any product in these publications.

For Research Use Only.

Not for use in diagnostic or therapeutic procedures. Not for resale. Not for distribution without written consent. Tonbo Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Tonbo Biosciences, Tonbo Biosciences Logo and all other trademarks are the property of Tonbo Biotechnologies Corporation. © 2013 Tonbo Biosciences.