

TECHNICAL DATA SHEET

APC Annexin V Apoptosis Detection Kit with 7-AAD

Catalog Number: 20-6410

PRODUCT INFORMATION

Kit Components APC Annexin V Cat.No. 20-6409-T100 Size 100 tests
 7-AAD Cat.No. 13-6993-T200 Size 200 tests
 Annexin V Binding Buffer (1X)

Description

Tonbo Biosciences APC Annexin V Apoptosis Detection Kit with 7-AAD contains the reagents needed identify apoptotic or necrotic cells by flow cytometry.

Tonbo Biosciences tests all antibodies by flow cytometry. Citations are provided as a resource for additional applications that have not been validated by Tonbo Biosciences. Please choose the appropriate format for each application and consult Materials and Methods sections for additional details about the use of any product in these publications.

For Research Use Only.

Not for use in diagnostic or therapeutic procedures. Not for resale. Not for distribution without written consent. Tonbo Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Tonbo Biosciences, Tonbo Biosciences Logo and all other trademarks are the property of Tonbo Biotechnologies Corporation. © 2013 Tonbo Biosciences.

TECHNICAL DATA SHEET

APC Annexin V

Catalog Number: 20-6409

PRODUCT INFORMATION


Contents: APC Annexin V

Concentration: 5µl/test

Use By: 12 months from date of receipt

Storage Conditions: 2-8°C protected from light

Formulation: 10 mM HEPES, 140 mM NaCl, 1mM EDTA, 0.09% NaN3 pH7.4

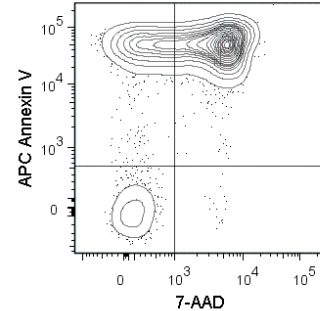

Mouse thymocytes were incubated overnight at 37°C in medium containing 1 µM dexamethasone. Cells were harvested and stained with APC Annexin V (20-6409) and 7-AAD Staining Solution (13-6993).

DESCRIPTION

Annexin V belongs to the annexin family of calcium-dependent phospholipid binding proteins and has a high affinity for phosphatidylserine (PS). One of the characteristic features of early stage apoptosis is the translocation of PS from the inner to the outer leaflet of the plasma membrane, exposing it to the extracellular environment where it can be detected. Fluorochrome conjugated Annexin V can serve as a sensitive detection probe for this hallmark of early stage apoptotic change.

Annexin V is generally used in combination with a viability dye, such as 7-AAD or propidium iodide (PI), that cannot pass through intact cellular membranes. Translocation of PS occurs prior to loss of membrane integrity, therefore viability dyes would be excluded from the cell. As cells progress through the apoptotic process, membranes are compromised and the viability dyes are able to enter the cell. Using this method, cells that are Annexin V positive and viability dye negative are in early apoptosis. If cells are positive for both Annexin V and the viability dye, they are either undergoing apoptotic death or necrosis. When negative for both markers, cells are viable.

PREPARATION & STORAGE

Annexin V was conjugated under optimal conditions, with unreacted dye removed from the preparation. It is recommended to store the product undiluted at 4°C, and protected from prolonged exposure to light. Do not freeze.

APPLICATION NOTES

This Annexin V preparation has been pre-titrated and quality-tested for flow cytometry using mouse thymocytes cultured overnight in the presence of dexamethasone to induce apoptosis. The Annexin V conjugate has been diluted for use at 5 µL per test, defined as the amount of protein that will stain a cell sample in a final volume of approximately 100 µL. The number of cells within a sample should be determined empirically, but typically ranges between 1x10⁵ to 1x10⁷ cells.

REFERENCES

Koopman G, Reutelingsperger CP, Kuijten GA, Keehnen RM, Pals ST, van Oers MH. 1994. *Blood*. 84(5):1415-1420. Vermes I, Haanen C, Steffens-Nakken H, Reutelingsperger C. 1995. *J. Immunol. Methods*. 184(1):39-51.

Tonbo Biosciences tests all antibodies by flow cytometry. Citations are provided as a resource for additional applications that have not been validated by Tonbo Biosciences. Please choose the appropriate format for each application and consult Materials and Methods sections for additional details about the use of any product in these publications.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale. Not for distribution without written consent. Tonbo Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Tonbo Biosciences, Tonbo Biosciences Logo and all other trademarks are the property of Tonbo Biotechnologies Corporation. © 2013 Tonbo Biosciences.

TECHNICAL DATA SHEET

7-AAD

Catalog Number: 13-6993

PRODUCT INFORMATION

Contents: 7-AAD

Excitation Laser Blue (488 nm), Green (532 nm), Yellow-Green (561 nm)

Emission (nm): 647

Formulation: 10 mM NaH₂PO₄, 150 mM NaCl, 0.09% NaN₃, pH 7.2

Storage Conditions: 2-8°C protected from light

Use by: 12 months from date of receipt

DESCRIPTION

7-AAD (7-Aminoactinomycin D) is a nucleic acid dye that can be used to exclude nonviable cells from the analysis of flow cytometry data.

PREPARATION & STORAGE

7-AAD is provided in solution and should be stored at 2-8°C and protected from light. Do not freeze. CAUTION: 7-AAD is a potential carcinogen. Protect skin and eyes by wearing suitable protective clothing, gloves and eye/face protection.

APPLICATION NOTES

7-AAD is useful as a viability probe for exclusion of nonviable cells based on light scatter and fluorescence properties. It is recommended to use 5 µL (0.25 µg) of 7-AAD solution per sample and incubate for 5-10 minutes in the dark prior to data acquisition.

REFERENCES

Schmid I, Krall WJ, Uittenbogaart CH, Braun J, Giorgi JV. Cytometry. 1992; 13(2):204-208.

NOTE: Please choose the appropriate format for each application. Citations are provided as a convenience to you; please consult Materials and Methods sections for additional details about the use of any product in these publications.

For Research Use Only.

Not for use in diagnostic or therapeutic procedures. Not for resale. Not for distribution without written consent. Tonbo Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Tonbo Biosciences, Tonbo Biosciences Logo and all other trademarks are the property of Tonbo Biotechnologies Corporation. © 2013 Tonbo Biosciences.

TECHNICAL DATA SHEET

Annexin V Binding Buffer (1X)

Catalog Number: TNB-5000

PRODUCT INFORMATION

Contents: Annexin V Binding Buffer (1X)

Use By: 6 months from date of receipt

Storage Conditions: 2-8°C

Formulation: Annexin V Binding Buffer (1X) contains HEPES, NaCl and CaCl₂.

DESCRIPTION

Annexin V Binding Buffer is formulated to facilitate the binding of Annexin V conjugates to phosphatidylserine for subsequent analysis by flow cytometry. This buffer is supplied as a 1X solution and requires no manipulation prior to use.

PREPARATION & STORAGE

Annexin V Binding Buffer (1X) is supplied as a 1X solution and requires no manipulation prior to use.

APPLICATION NOTES

Annexin V Binding Buffer (1X) is recommended for use with Annexin V reagents. Cells should be resuspended in Annexin V Binding Buffer (1X) prior to incubation with Annexin V reagents.

Tonbo Biosciences tests all antibodies by flow cytometry. Citations are provided as a resource for additional applications that have not been validated by Tonbo Biosciences. Please choose the appropriate format for each application and consult Materials and Methods sections for additional details about the use of any product in these publications.

For Research Use Only.

Not for use in diagnostic or therapeutic procedures. Not for resale. Not for distribution without written consent. Tonbo Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Tonbo Biosciences, Tonbo Biosciences Logo and all other trademarks are the property of Tonbo Biotechnologies Corporation. © 2013 Tonbo Biosciences.