

## TECHNICAL DATA SHEET

# APC Anti-Mouse IL-12/IL-23 p40 (C17.8)

Catalog Number: 20-7123

## PRODUCT INFORMATION

**Contents:** APC Anti-Mouse IL-12/IL-23 p40 (C17.8)

**Isotype:** Rat IgG2a, kappa

**Concentration:** 0.2 mg/mL


**Clone:** C17.8

**Reactivity:** Mouse

**Use By:** 12 months from date of receipt

**Storage Conditions:** 2-8°C protected from light

**Formulation:** 10 mM NaH<sub>2</sub>PO<sub>4</sub>, 150 mM NaCl, 0.09% NaN<sub>3</sub>,  
0.1% gelatin, pH7.2

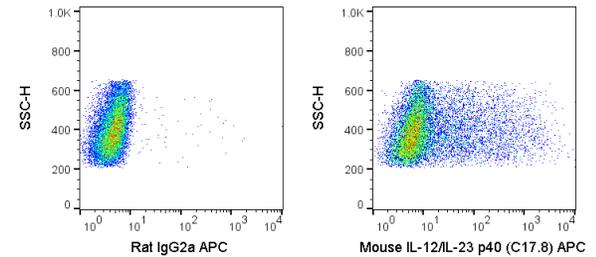

Mouse macrophages were stimulated in the presence of a protein transport inhibitor. Cells were then fixed, permeabilized and stained intracellularly with 0.06 ug APC Anti-Mouse IL-12/IL-23 p40 (20-7123) (left panel) or 0.06 ug APC Rat IgG2a (left panel).

## DESCRIPTION

The C17.8 antibody is specific for the 40 kDa (p40) protein subunit shared by the cytokines IL-12 and IL-23. To form IL-12, p40 assembles with a separate 35 kDa protein known as p35, resulting in a 70 kDa functional cytokine. IL-12 is secreted by activated monocytes, macrophages, and dendritic cells, and has been shown to target naive, resting CD4+ T cells to promote their proliferation and secretion of cytokines. IL-23 contains the p40 subunit in combination with a 19 kDa protein chain, p19; its primary source being activated dendritic cells and other antigen-presenting cells. IL-23 appears to target different cell types than IL-12, acting on memory CD4+ T cells to induce a strong proliferative response and contributing to the generation and expansion of Th17 cells. Like the cytokines themselves, the receptors for IL-12 and IL-23 share one subunit, as well as containing distinct cytokine-specific subunits. As the C17.8 antibody binds to a shared subunit of both IL-12 and IL-23, it may be used as a marker for either IL-12 or IL-23 expression in dendritic cells, monocytes and macrophages, and is widely used for neutralization of activity associated with either cytokine. Please choose the appropriate format for each application.

## PREPARATION & STORAGE

This monoclonal antibody was purified from tissue culture supernatant via affinity chromatography. The purified antibody was conjugated under optimal conditions, with unreacted dye removed from the preparation. It is recommended to store the product undiluted at 4°C, and protected from prolonged exposure to light. Do not freeze.

## APPLICATION NOTES

This antibody preparation has been quality-tested for flow cytometry using mouse spleen cells, or an appropriate cell type (where indicated). Please refer to the figure legend for the optimal concentration used to stain the tissue shown. We recommend titrating the antibody under your specific conditions to determine the optimal concentration of antibody needed in your experimental system.

## REFERENCES

Dong H, Franklin NA, Roberts DJ, Yagita H, Glennie MJ and Bullock TNJ. 2012. J. Immunol. 188: 3829-3838. (in vivo blocking) Prabhakara R, Harro JM, Leid JG, Keegan AD, Prior ML, and Shirliff ME. 2011. Infect. Immun. 79: 5010-5018. (in vivo blocking) Chmielewski M, Kopecky C, Hombach AA, and Abken H. 2011. Cancer Res. 71: 5697-5706. (ELISA) Lo C-H, Lee S-C, Wu P-Y, Pan W-Y et al. 2003. J. Immunol. 171: 600-607. (Immunoprecipitation) Belladonna ML, Renaud J-C, Bianchi R, Vacca C, Fallarino F, Orabona C, Fioretti MC, Grohmann, and Puccetti P. 2002. J. Immunol. 168: 5448-5454. (Western Blot) Ludviksson BR, Ehrhardt RO, and Strober W. 1999. J. Immunol. 163: 4349-4359. (Immunofluorescence microscopy – frozen tissue) Kato K, Shimozato O, Hoshi K, Wakimoto H, Hamada H, Yagita H, and Okumura K. 1996. Proc. Natl. Acad. Sci. 93:9085-9089. (Immunoprecipitation; ELISA) Wysocka M, Kubin M, Vieira LQ, Ozmen L, Garotta G, Scott P and Trinchieri G. 1995. Eur. J. Immunol. 25: 672-676. (Origination of clone; Immunoprecipitation, Western Blot)

Tonbo Biosciences tests all antibodies by flow cytometry. Citations are provided as a resource for additional applications that have not been validated by Tonbo Biosciences. Please choose the appropriate format for each application and consult Materials and Methods sections for additional details about the use of any product in these publications.

For Research Use Only.

Not for use in diagnostic or therapeutic procedures. Not for resale. Not for distribution without written consent. Tonbo Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Tonbo Biosciences, Tonbo Biosciences Logo and all other trademarks are the property of Tonbo Biotechnologies Corporation. © 2013 Tonbo Biosciences.