

TECHNICAL DATA SHEET

Biotin Anti-Human CD56 (NCAM) (MY31)

Catalog Number: 30-0564

PRODUCT INFORMATION

Contents: Biotin Anti-Human CD56 (NCAM) (MY31)

Isotype: Mouse IgG1, kappa

Concentration: 0.5 mg/mL


Clone: MY31

Reactivity: Human

Use By: 12 months from date of receipt

Storage Conditions: 2-8°C

Formulation: 10 mM NaH₂PO₄, 150 mM NaCl, 0.09% NaN₃, pH 7.2

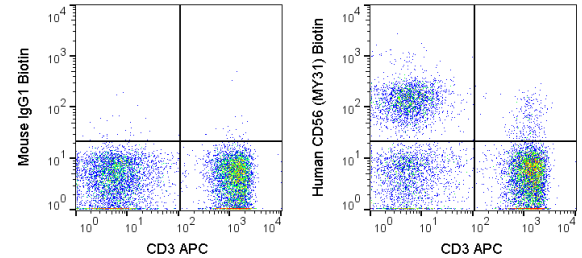

Human peripheral blood lymphocytes were stained with APC Anti-Human CD3 (20-0038) and 0.5 ug Biotin Anti-Human CD56 (30-0564) (right panel) or 0.5 ug Biotin Mouse IgG1 isotype control (left panel), followed by Streptavidin PE.

DESCRIPTION

The MY31 antibody reacts with human CD56, also known as the neural cell adhesion molecule (NCAM), a glycoprotein which is a member of the immunoglobulin superfamily. The 140 kDa isoform of CD56 is expressed on human NK cells and NK-T cells, with increased expression levels on activated NK lymphocytes. The CD56 antigen is also expressed by neurons and is reported to play a role in nervous system development and neural cell-to-cell adhesion. Clone MY31 also reacts with a subset of CD14+ monocytes in non-human primates, and is reported to be cross-reactive with Chimpanzee, Cynomolgus and Rhesus.

PREPARATION & STORAGE

This monoclonal antibody was purified from tissue culture supernatant via affinity chromatography. The purified antibody was conjugated under optimal conditions, with unreacted biotin removed from the preparation. It is recommended to store the product undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

APPLICATION NOTES

This antibody preparation has been quality-tested for flow cytometry using an appropriate cell type (as indicated). Please refer to the figure legend for the optimal concentration used to stain the tissue shown. We recommend titrating the antibody under your specific conditions to determine the optimal concentration of antibody needed in your experimental system.

REFERENCES

- Lanier LL, Testi R, Bindl J and Phillips JH. 1989. *J Exp Med.* 169: 2233-2238.
- Schlossman SF, Boumsell L, Gilks W et al., eds. 1995. *Leucocyte Typing V: White Cell Differentiation Antigens.* Oxford University Press.
- Carter DL, Shieh TM, Blosser RL, Chadwick KR, Margolick JB, Hidreth JE, Clements JE and Zink MC. 1999. *Cytometry.* 37(1): 41-50.
- Chan WK, Suwannasaen D, Throm RE, Li Y, Eldridge PW, Houston J, Gray JT, Pui C-H and Leung W. 2015. *Leukemia.* 29: 387-395. (Flow cytometry)
- Woltman AM, Op den Brouw ML, Biesta PJ, Shi CC and Janssen HLA. 2011. *PLoS ONE* 6(1): e15324. doi: 10.1371/journal.pone.0015324. (Flow cytometry)
- Bleul CC, Wu L, Hoxie JA, Springer TA and Mackay CR. 1997. *Proc Natl Acad Sci USA.* 94(5): 1925-1930. (Flow cytometry)
- Brown K and Barratt-Boyes SM. 2009. *J Med Primatol.* 38(4): 272-278. (Flow cytometry – Rhesus)

Tonbo Biosciences tests all antibodies by flow cytometry. Citations are provided as a resource for additional applications that have not been validated by Tonbo Biosciences. Please choose the appropriate format for each application and consult Materials and Methods sections for additional details about the use of any product in these publications.

For Research Use Only.

Not for use in diagnostic or therapeutic procedures. Not for resale. Not for distribution without written consent. Tonbo Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Tonbo Biosciences, Tonbo Biosciences Logo and all other trademarks are the property of Tonbo Biotechnologies Corporation. © 2013 Tonbo Biosciences.