

TECHNICAL DATA SHEET

# FITC Anti-Mouse CD5 (53-7.3)

Catalog Number: 35-0051

## PRODUCT INFORMATION

**Contents:** FITC Anti-Mouse CD5 (53-7.3)

**Isotype:** Rat IgG2a, kappa

**Concentration:** 0.5 mg/mL


**Clone:** 53-7.3

**Reactivity:** Mouse

**Use By:** 12 months from date of receipt

**Storage Conditions:** 2-8°C protected from light

**Formulation:** 10 mM NaH<sub>2</sub>PO<sub>4</sub>, 150 mM NaCl, 0.09% NaN<sub>3</sub>, 0.1% gelatin, pH7.2

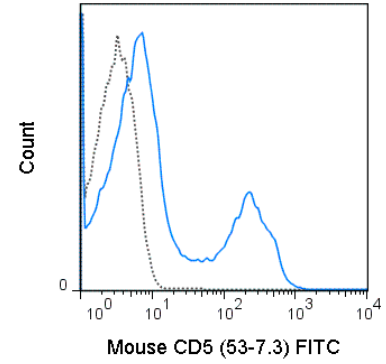

C57Bl/6 splenocytes were stained with 0.5 ug FITC Anti-Mouse CD5 (35-0051) (solid line) or 0.5 ug FITC Rat IgG2a isotype control (dashed line).

## DESCRIPTION

The 53-7.3 monoclonal antibody is specific for mouse CD5, a 67 kDa protein that is expressed on most thymocytes, mature T cells and the B-1 subset of B cells. CD5 is a member of the scavenger receptor cysteine-rich protein superfamily and is a ligand for the c-type lectin CD72. CD5 is also known as Leu-1 and functions to modulate antigen receptor signaling on T and B cells.

## PREPARATION & STORAGE

This monoclonal antibody was purified from tissue culture supernatant via affinity chromatography. The purified antibody was conjugated under optimal conditions, with unreacted dye removed from the preparation. It is recommended to store the product undiluted at 4°C, and protected from prolonged exposure to light. Do not freeze.

## APPLICATION NOTES

This antibody preparation has been quality-tested for flow cytometry using mouse spleen cells, or an appropriate cell type (where indicated). Please refer to the figure legend for the optimal concentration used to stain the tissue shown. We recommend titrating the antibody under your specific conditions to determine the optimal concentration of antibody needed in your experimental system.

## REFERENCES

Azzam HS, Grinberg A, Lui K, Shen H, Shores EW, Love PE. 1998. J Exp Med. 188(12):2301-2311. Bikah G, Carey J, Ciallella JR, Tarakhovsky A, Bondada S. 1996. Science. 274(5294):1906-1909. Bikah G, Lynd FM, Aruffo AA, Ledbetter JA, Bondada S. 1998. Int Immunol. 10(8):1185-1196. Hayakawa K, Hardy RR, Parks DR, Herzenberg LA. 1983. J Exp Med. 157(1):202-218. Ledbetter JA, Rouse RV, Micklem HS, Herzenberg LA. 1980. J Exp Med. 152(2):280-295. Tarakhovsky A, Kanner SB, Hombach J, et al. 1995. Science. 269(5223):535-537.

Tonbo Biosciences tests all antibodies by flow cytometry. Citations are provided as a resource for additional applications that have not been validated by Tonbo Biosciences. Please choose the appropriate format for each application and consult Materials and Methods sections for additional details about the use of any product in these publications.

For Research Use Only.

Not for use in diagnostic or therapeutic procedures. Not for resale. Not for distribution without written consent. Tonbo Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Tonbo Biosciences, Tonbo Biosciences Logo and all other trademarks are the property of Tonbo Biotechnologies Corporation. © 2013 Tonbo Biosciences.