

TECHNICAL DATA SHEET

PE Anti-Human CD27 (O323)

Catalog Number: 50-0279

PRODUCT INFORMATION

Contents: PE Anti-Human CD27 (O323)


Isotype: Mouse IgG1, kappa

Concentration: 5 uL (0.25 ug)/test

Clone: O323

Reactivity: Human

Formulation: 10 mM NaH₂PO₄, 150 mM NaCl, 0.09% NaN₃,
0.1% gelatin, pH7.2

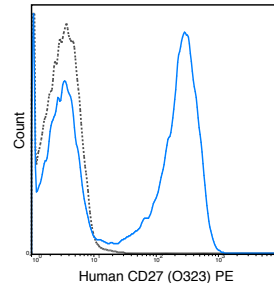

Human peripheral blood lymphocytes were stained with 5 uL (0.25 ug) PE Anti-Human CD27 (50-0279) (solid line) or 0.25 ug PE Mouse IgG1 isotype control.

DESCRIPTION

The O323 antibody reacts with human CD27 (TNFRSF7), a cell surface homodimer of 55 kDa subunits, which provides co-stimulatory signaling in support of the T cell (TCR) and B cell (BCR) receptors. By comparison with CD28, whose TCR co-stimulatory signal can trigger cell proliferation, CD27 signaling appears to promote cell survival and differentiation to effector / memory stages. Also in contrast with CD28, the CD27 receptor may be shed following interaction with its ligand CD70, which is typically expressed on activated dendritic cells, T cells and B cells. With respect to B cells, CD27 is considered to be a phenotypic marker for memory B cells. CD27 has been included within a group of phenotypic markers for identifying human B regulatory cells (Bregs), a cell type proposed to regulate CD4+ T cell proliferation and Foxp3 / CTLA-4 expression in Treg cells. The O323 antibody may be used for analysis of CD27 expression on peripheral T cells, and is frequently used in combination with antibodies for IgD and IgM to distinguish naïve vs. memory B cell populations. For identification of Breg cells, this antibody has been used in combination with antibodies for CD25, CD27, CD1d, IL-10 and TGF-beta (Kessel et al. 2012. *Autoimm. Rev.* 11(9): 670-677). The antibody is also reported for cross-reactivity with Baboon, Cynomolgus and Rhesus CD27.

PREPARATION & STORAGE

This monoclonal antibody was purified from tissue culture supernatant via affinity chromatography. The purified antibody was conjugated under optimal conditions, with unreacted dye removed from the preparation. It is recommended to store the product undiluted at 4°C, and protected from prolonged exposure to light. Do not freeze.

APPLICATION NOTES

This antibody preparation has been pre-titrated and quality-tested for flow cytometry using an appropriate cell type. The antibody has been diluted for use at 5 uL per test, defined as the amount of antibody that will stain a cell sample in a final volume of approximately 100 uL. The number of cells within a sample should be determined empirically, but typically ranges between 1x10⁵ to 1x10⁸ cells.

REFERENCES

Kroenke MA, Eto D, Locci M, Cho M, Davidson T, Haddad EK, and Crotty S. 2012. *J. Immunol.* 188: 3734-3744. (flow cytometry) Ruffell B, Au A, Rugo HS, Esserman LJ, Hwang ES, and Coussens LM. 2012. *Proc. Natl. Acad. Sci.* 109: 2796-2801. (flow cytometry) Strowig T, Tongvaux A, Rathinam C, Takizawa H, Borsotti C, Philbrick W, Eynon EE, Manz MG, and Flavell RA. 2011. *Proc. Natl. Acad. Sci.* 108: 13218-13223. (flow cytometry) So NSY, Ostrowski MA and Gray-Owen SD. 2012. *J. Immunol.* 188: 4008-4022. (flow cytometry) Klatt NR, Vinton CL, Lynch RM, Canary LA, Ho J, Darrah PA, Estes JD, Seder RA, Moir SL, and Brenchley JM. 2011. *Blood.* 118: 5803-5812. (Cell sorting – Rhesus) Ali Z, Yan L, Plagman N, Reichenberg A, Hintz M, Jomaa H, Villinger F, and Chen ZW. 2009. *J. Immunol.* 183: 5407-5417. (flow cytometry – Rhesus) Ali Z, Shao L, Halliday L, Reichenberg A, Hintz M, Jomaa H, and Chen ZW. 2007. *J. Immunol.* 8287-8296. (flow cytometry – Cynomolgus)

NOTE: Please choose the appropriate format for each application. Citations are provided as a convenience to you; please consult Materials and Methods sections for additional details about the use of any product in these publications.

For Research Use Only.

Not for use in diagnostic or therapeutic procedures. Not for resale. Not for distribution without written consent. Tonbo Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Tonbo Biosciences, Tonbo Biosciences Logo and all other trademarks are the property of Tonbo Biotechnologies Corporation. © 2013 Tonbo Biosciences.