

TECHNICAL DATA SHEET

# PE-Cy7 Anti-Human CD19 (SJ25C1)

Catalog Number: 60-0198

## PRODUCT INFORMATION

**Contents:** PE-Cy7 Anti-Human CD19 (SJ25C1)

**Isotype:** Mouse IgG1, kappa

**Concentration:** 5 µL (0.125 µg)/test


**Clone:** SJ25C1

**Reactivity:** Human

**Use By:** 6 months from date of receipt

**Storage Conditions:** 2-8°C protected from light

**Formulation:** 10 mM NaH<sub>2</sub>PO<sub>4</sub>, 150 mM NaCl, 0.09% NaN<sub>3</sub>, 0.1% gelatin, pH7.2

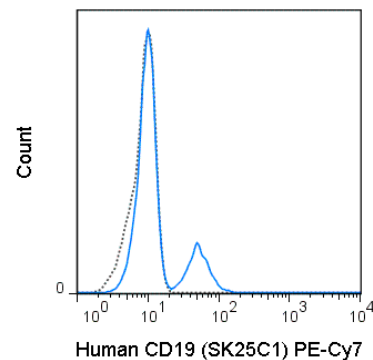

Human peripheral blood lymphocytes were stained with 5 µL (0.125 µg) PE-Cy7 Anti-Human CD19 (60-0198) (solid line) or 0.125 µg PE-Cy7 Mouse IgG1 isotype control (dashed line).

## DESCRIPTION

The SJ25C1 antibody reacts with human CD19, a 95 kDa glycoprotein which acts as a co-receptor, along with CD21 (CR2), CD81 (TAPA-1) and CD225 (Leu13), in support of the functional B cell receptor (BCR). This complex provides antigen-specific recognition and subsequent activation of B cells to proliferate and differentiate into antibody-secreting cells (plasma cells) or memory B cells, which are crucial for secondary antigen encounter. Upon activation and tyrosine phosphorylation, the CD19 molecule can provide an anchor for cytoplasmic signaling proteins such as GRB2, SOS or PLCG2. CD19 is a lineage-differentiation marker, as its expression is detectable at the earliest B cell stages, through development, and is finally lost upon transition to mature plasma cells. The SJ25C1 antibody is widely used as a phenotypic marker for CD19 expression on B cells, as well as on dendritic cell subsets.

## PREPARATION & STORAGE

This monoclonal antibody was purified from tissue culture supernatant via affinity chromatography. The purified antibody was conjugated under optimal conditions, with unreacted dye removed from the preparation. It is recommended to store the product undiluted at 4°C, and protected from prolonged exposure to light. Do not freeze.

## APPLICATION NOTES

This antibody preparation has been pre-titrated and quality-tested for flow cytometry using an appropriate cell type. The antibody has been diluted for use at 5 µL per test, defined as the amount of antibody that will stain a cell sample in a final volume of approximately 100 µL. The number of cells within a sample should be determined empirically, but typically ranges between 1x10<sup>5</sup> to 1x10<sup>8</sup> cells.

## REFERENCES

- Piatosa B, Birbach M, Siewiera K, Ussowicz M, Kalwak K, Drabko K, Rekawek A, Tkaczyk K, Kurowski PN. 2013. *Cytometry Part B*. 84B: 179–186. (Flow Cytometry)
- Sekiguchi DR, Smith SB, Sutter JA, Goodman NG, Probert K, Louzoun Y, Rogers W, Luning Prak ET. 2011. *Cytometry Part B*. 80B: 291–299. (Flow Cytometry)
- Dörken B, Möller P, Pezzutto A, Schwartz-Albiez R, Moldenhauer G. B-cell antigens: CD19. In: Knapp W, Dörken B, Gilks WR, et al, ed. *Leucocyte Typing IV: White Cell Differentiation Antigens*. New York, NY: Oxford University Press; 1989:34-36. (Flow Cytometry)
- Karahan GE, Eikmans M, Anholts JDH, Class FHJ and Heidt S. 2014. *Clin Exp Immunol*. 177(1): 333-340. (Flow Cytometry)
- Bansal RR, Mackay CR, Moser B and Ebert M. 2012. *Eur J Immunol*. 42(1): 110-119. (Flow Cytometry)
- Jourdan M, Caraux A, Caron G, Robert N, Fiol G, Rème T, Bolloré K, Vendrell JP, Le Gallou S, Mourcin F, De Vos J, Kassambara A, Duperray C, Hose D, Fest T, Tarte K and Klein B. 2011. *J Immunol*. 187(8): 3931-3941. (Flow Cytometry)
- Goval J-J, Greimers R, Boniver J and de Leval L. 2006. *J Histochem Cytochem*. 54(1): 75-84. (Immunofluorescence – frozen tissue)

Tonbo Biosciences tests all antibodies by flow cytometry. Citations are provided as a resource for additional applications that have not been validated by Tonbo Biosciences. Please choose the appropriate format for each application and consult Materials and Methods sections for additional details about the use of any product in these publications.

For Research Use Only.

Not for use in diagnostic or therapeutic procedures. Not for resale. Not for distribution without written consent. Tonbo Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Tonbo Biosciences, Tonbo Biosciences Logo and all other trademarks are the property of Tonbo Biotechnologies Corporation. © 2013 Tonbo Biosciences.