

TECHNICAL DATA SHEET

PerCP-Cy5.5 Anti-Human CD5 (L17F12)

Catalog Number: 65-0058

PRODUCT INFORMATION

Contents: PerCP-Cy5.5 Anti-Human CD5 (L17F12)

Isotype: Mouse IgG2a, kappa

Concentration: 5 µL (0.06 µg)/test


Clone: L17F12

Reactivity: Human

Use By: 6 months from date of receipt

Storage Conditions: 2-8°C protected from light

Formulation: 10 mM NaH₂PO₄, 150 mM NaCl, 0.09% NaN₃, 0.1% gelatin, pH7.2

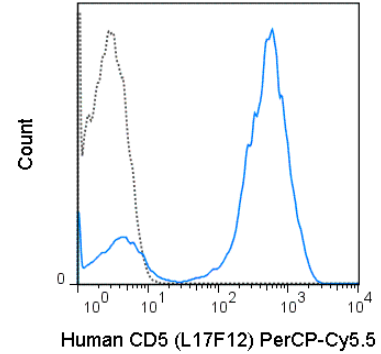

Human peripheral blood lymphocytes were stained with 5 µL (0.06 µg) PerCP-Cy5.5 Anti-Human CD5 (65-0058) (solid line) or 0.06 µg PerCP-Cy5.5 Mouse IgG2a isotype control (dashed line).

DESCRIPTION

The L17F12 antibody is specific for human CD5, a 67 kD transmembrane glycoprotein that is expressed on most thymocytes, mature T cells, and a subset of B cells. CD5 is a member of the scavenger receptor superfamily and is present on approximately 70% of normal peripheral blood lymphocytes. CD5 is involved in modulating antigen receptor signaling in both T and B cells and acts through binding of CD72, its receptor expressed on B cells.

PREPARATION & STORAGE

This monoclonal antibody was purified from tissue culture supernatant via affinity chromatography. The purified antibody was conjugated under optimal conditions, with unreacted dye removed from the preparation. It is recommended to store the product undiluted at 4°C, and protected from prolonged exposure to light. Do not freeze.

APPLICATION NOTES

This antibody preparation has been pre-titrated and quality-tested for flow cytometry using an appropriate cell type. The antibody has been diluted for use at 5 µL per test, defined as the amount of antibody that will stain a cell sample in a final volume of approximately 100 µL. The number of cells within a sample should be determined empirically, but typically ranges between 1x10⁵ to 1x10⁸ cells.

REFERENCES

Knapp W, Dorken B, et al. eds. 1989. Leucocyte Typing IV: White Cell Differentiation Antigens. Oxford University Press. New York.
 Gadol N and Ault KA. 1986. Immunol Rev. 93: 23-34.
 Antin JH, Emerson SG, Martin P, Gadol N and Ault KA. 1986. J Immunol. 136: 505-510.
 Raab M, Yamamoto M and Rudd CE. 1994. Mol Cell Biol. 14(5): 2862-2870.
 Lydyard PM, Lamour A, MacKenzie LE, Jamin C, Magged RA and Youinou P. 1993. Immunol Lett. 38(2): 159-166.
 McKenna RW, Washington LT, Aquino DB, Picker LJ and Kroft SH. 2001. Blood. 98(8): 2498-2507. (Flow Cytometry)
 Butmarc JR, Kourea HP, Levi E and Kadin ME. 1998. Am J Clin Path. 109(6): 682-688. (Immunohistochemistry - paraffin-embedded tissue)

Tonbo Biosciences tests all antibodies by flow cytometry. Citations are provided as a resource for additional applications that have not been validated by Tonbo Biosciences. Please choose the appropriate format for each application and consult Materials and Methods sections for additional details about the use of any product in these publications.

For Research Use Only.

Not for use in diagnostic or therapeutic procedures. Not for resale. Not for distribution without written consent. Tonbo Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Tonbo Biosciences, Tonbo Biosciences Logo and all other trademarks are the property of Tonbo Biotechnologies Corporation. © 2013 Tonbo Biosciences.