

TECHNICAL DATA SHEET

# PerCP-Cyanine5.5 Anti-Human CD38 (HIT2)

Catalog Number: 65-0389

## PRODUCT INFORMATION

**Contents:** PerCP-Cyanine5.5 Anti-Human CD38 (HIT2)

**Isotype:** Mouse IgG1, kappa

**Concentration:** 5ul (0.25 ug)/test


**Clone:** HIT2

**Reactivity:** Human

**Use By:** 6 months from date of receipt

**Storage Conditions:** 2-8°C protected from light

**Formulation:** 10 mM NaH<sub>2</sub>PO<sub>4</sub>, 150 mM NaCl, 0.09% NaN<sub>3</sub>, 0.1% gelatin, pH7.2

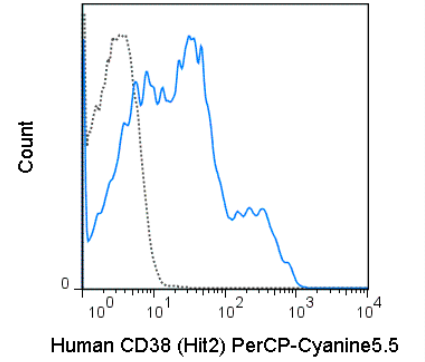

Human peripheral blood lymphocytes were stained with 5 uL (0.25 ug) PerCP-Cyanine5.5 Anti-Human CD38 (65-0389) (solid line) or 0.25 ug PerCP-Cyanine5.5 Mouse IgG1 isotype control (dashed line).

## DESCRIPTION

The HIT2 antibody is specific for human CD38, a 45 kD type II transmembrane glycoprotein expressed on thymocytes, activated T and B cells, and monocytes as well as other non-hematopoietic cells. CD38 is an ectoenzyme that functions to catalyze the synthesis and hydrolysis of cyclic ADP-ribose and is involved in cell signaling, regulating cell adhesion, and activation. Additionally, CD38 has been shown to be a prognostic marker for some leukemias and other diseases.

## PREPARATION & STORAGE

This monoclonal antibody was purified from tissue culture supernatant via affinity chromatography. The purified antibody was conjugated under optimal conditions, with unreacted dye removed from the preparation. It is recommended to store the product undiluted at 4°C, and protected from prolonged exposure to light. Do not freeze.

## APPLICATION NOTES

This antibody preparation has been pre-titrated and quality-tested for flow cytometry using an appropriate cell type. The antibody has been diluted for use at 5 µL per test, defined as the amount of antibody that will stain a cell sample in a final volume of approximately 100 µL. The number of cells within a sample should be determined empirically, but typically ranges between 1x10e5 to 1x10e8 cells.

## REFERENCES

- Dianzani U, and Malavasi F. 1995. Crit. Rev. Immunol. 15(2):167-200.
- Mehta K, Shahid U, and Malavasi F. Oct. 1996. FASEB J. 10(12):1408-1417.
- Ferrero E, and Malavas F. Feb. 1999. J. Leukoc. Biol. 65(2):151-161.
- Matrai Z. Feb. 2005. Hematology. 10(1):39-46.

Tonbo Biosciences tests all antibodies by flow cytometry. Citations are provided as a resource for additional applications that have not been validated by Tonbo Biosciences. Please choose the appropriate format for each application and consult Materials and Methods sections for additional details about the use of any product in these publications.

For Research Use Only.

Not for use in diagnostic or therapeutic procedures. Not for resale. Not for distribution without written consent. Tonbo Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Tonbo Biosciences, Tonbo Biosciences Logo and all other trademarks are the property of Tonbo Biotechnologies Corporation. © 2013 Tonbo Biosciences.