

TECHNICAL DATA SHEET

PerCP-Cyanine5.5 Anti-Human CD69 (FN50)

Catalog Number: 65-0699

PRODUCT INFORMATION

Contents: PerCP-Cyanine5.5 Anti-Human CD69 (FN50)

Isotype: Mouse IgG1, k

Concentration: 5 uL (0.5 ug)/test


Clone: FN50

Reactivity: Human

Use By: 6 months from date of receipt

Storage Conditions: 2-8°C protected from light

Formulation: 10 mM NaH₂PO₄, 150 mM NaCl, 0.09% NaN₃, 0.1% gelatin, pH7.2

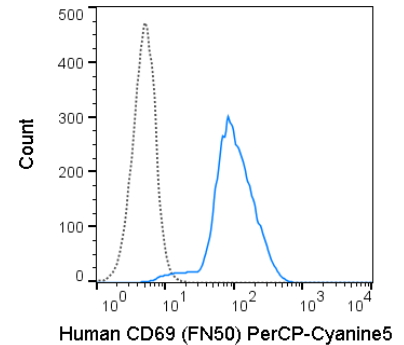

Human peripheral blood lymphocytes were stimulated with PMA+Ionomycin for 5 hours and then stained with 5 uL (0.5 ug) PerCP-Cyanine5.5 Anti-Human CD69 (65-0699) (solid line) or 0.5 ug PerCP-Cyanine5.5 Mouse IgG1 isotype control (dashed line).

DESCRIPTION

The FN50 monoclonal antibody reacts with human CD69, a type II transmembrane glycoprotein also known as the Very Early Activation Antigen, EA-1, Leu23, Activation Inducer Molecule (AIM) and CLEC2C. CD69 is expressed as a disulfide-linked homodimer on activated T and B cells, NK cells, neutrophils, and monocytes. Induction occurs rapidly upon activation and is transient. It is also constitutively expressed on platelets and a subset of thymocytes.

PREPARATION & STORAGE

This monoclonal antibody was purified from tissue culture supernatant via affinity chromatography. The purified antibody was conjugated under optimal conditions, with unreacted dye removed from the preparation. It is recommended to store the product undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

APPLICATION NOTES

This antibody preparation has been pre-titrated and quality-tested for flow cytometry using an appropriate cell type. The antibody has been diluted for use at 5 uL per test, defined as the amount of antibody that will stain a cell sample in a final volume of approximately 100 uL. The number of cells within a sample should be determined empirically, but typically ranges between 1x10e5 to 1x10e8 cells.

REFERENCES

Testi R, D'Ambrosio D, De Maria R, Santoni A. 1994. Immunol Today. Oct;15(10):479-483. Schlossman S, et al. Eds. 1995. Leucocyte Typing V. Oxford University Press. New York.

Tonbo Biosciences tests all antibodies by flow cytometry. Citations are provided as a resource for additional applications that have not been validated by Tonbo Biosciences. Please choose the appropriate format for each application and consult Materials and Methods sections for additional details about the use of any product in these publications.

For Research Use Only.

Not for use in diagnostic or therapeutic procedures. Not for resale. Not for distribution without written consent. Tonbo Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Tonbo Biosciences, Tonbo Biosciences Logo and all other trademarks are the property of Tonbo Biotechnologies Corporation. © 2013 Tonbo Biosciences.