

TECHNICAL DATA SHEET

violetFluor™ 450 Anti-Human CD38 (HIT2)

Catalog Number: 75-0389

PRODUCT INFORMATION

Contents: violetFluor™ 450 Anti-Human CD38

Isotype: Mouse IgG1, kappa

Concentration: 5 uL (0.5 ug)/test


Clone: HIT2

Reactivity: Human

Use By: 12 months from date of receipt

Storage Conditions: 2-8°C protected from light

Formulation: 10 mM NaH₂PO₄, 150 mM NaCl,
0.09% NaN₃, 0.1% gelatin, pH7.2

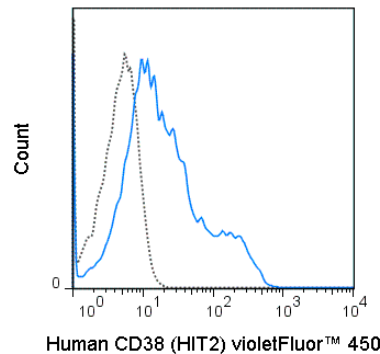

Human peripheral blood lymphocytes were stained with 5 uL (0.5 ug) violetFluor™ 450 Anti-Human CD38 (75-0389) (solid line) or 0.5 ug violetFluor™ 450 Mouse IgG1 isotype control (dashed line).

DESCRIPTION

The HIT2 antibody is specific for human CD38, a 45 kD type II transmembrane glycoprotein expressed on thymocytes, activated T and B cells, and monocytes as well as other non-hematopoietic cells. CD38 is an ectoenzyme that functions to catalyze the synthesis and hydrolysis of cyclic ADP-ribose and is involved in cell signaling, regulating cell adhesion, and activation. Additionally, CD38 has been shown to be a prognostic marker for some leukemias and other diseases.

PREPARATION & STORAGE

This monoclonal antibody was purified from tissue culture supernatant via affinity chromatography. The purified antibody was conjugated under optimal conditions, with unreacted dye removed from the preparation. It is recommended to store the product undiluted at 4°C, and protected from prolonged exposure to light. Do not freeze.

APPLICATION NOTES

This antibody preparation has been pre-titrated and quality-tested for flow cytometry using an appropriate cell type. The antibody has been diluted for use at 5 uL per test, defined as the amount of antibody that will stain a cell sample in a final volume of approximately 100 uL. The number of cells within a sample should be determined empirically, but typically ranges between 1x10⁵ to 1x10⁸ cells.

REFERENCES

Dianzani U, and Malavasi F. 1995. Crit. Rev. Immunol. 15(2):167-200. Mehta K, Shahid U, and Malavasi F. Oct. 1996. FASEB J. 10(12):1408-1417. Ferrero E, and Malavas F. Feb. 1999. J. Leukoc. Biol. 65(2):151-161. Matrai Z. Feb. 2005. Hematology. 10(1):39-46.

NOTE: Please choose the appropriate format for each application. Citations are provided as a convenience to you; please consult Materials and Methods sections for additional details about the use of any product in these publications.

For Research Use Only.

Not for use in diagnostic or therapeutic procedures. Not for resale. Not for distribution without written consent. Tonbo Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Tonbo Biosciences, Tonbo Biosciences Logo and all other trademarks are the property of Tonbo Biotechnologies Corporation. © 2013 Tonbo Biosciences.