

TECHNICAL DATA SHEET

violetFluor™ 450 Anti-Human/Mouse CD45R (B220) (RA3-6B2)

Catalog Number: 75-0452

PRODUCT INFORMATION

Contents: violetFluor™ 450 Anti-Human/Mouse CD45R (B220) (RA3-6B2)


Isotype: Rat IgG2a, kappa

Concentration: 1 mg/ml

Clone: RA3-6B2

Reactivity: Human, Mouse

Formulation: 10 mM NaH₂PO₄, 150 mM NaCl, 0.09% NaN₃, 0.1% gelatin, pH7.2

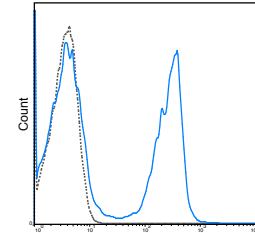

Hu/Mo CD45R (B220) (RA3-6B2) violetFluor 450

C57Bl/6 splenocytes were stained with 0.5 ug violetFluor 450 Anti-Hu/Mo CD45R (B220) (75-0452) (solid line) or 0.5 ug violetFluor 450 Rat IgG2a isotype control (dashed line).

DESCRIPTION

The RA3-6B2 antibody reacts with the human and mouse CD45 isoform known as CD45R, or B220, a protein tyrosine phosphatase of 220 kDa. CD45 is one of the most abundant hematopoietic markers, and is expressed on all leukocytes (the Leukocyte Common Antigen, LCA). Various isoforms are generated and expressed in cell-specific patterns, all critical for leukocyte function. In mouse, the CD45R/B220 isoform is predominantly found on B cells, at varying levels on all stages from pro-B cells to activated B cells, and may also be detected on certain T cell and NK cell subsets. It is of note that B220 is not similarly expressed on human B cells, where it appears to be differentiation-specific and therefore expressed on only some B cell subsets. Other forms of CD45 with restricted cellular expression include CD45RA, CD45RB, CD45RO and several others. The RA3-6B2 antibody is one of the most consistently used leukocyte markers for B cells, T cell subsets and NK cell subsets in human and mouse. This antibody has also been reported as cross-reactive with feline CD45R/B220.

PREPARATION & STORAGE

This monoclonal antibody was purified from tissue culture supernatant via affinity chromatography. The purified antibody was conjugated under optimal conditions, with unreacted dye removed from the preparation. It is recommended to store the product undiluted at 4°C, and protected from prolonged exposure to light. Do not freeze.

APPLICATION NOTES

This antibody preparation has been quality-tested for flow cytometry using mouse spleen cells, or an appropriate cell type (where indicated). The amount of antibody required for optimal staining of a cell sample should be determined empirically in your system.

violetFluor™ 450 dye is excited by the violet (405 nm) laser and has a peak emission of 450 nm. The most common band pass filters for this dye are 440/40 or 450/50. violetFluor™ 450 can be used as an alternative for Pacific Blue®, BD Horizon™ V450 or eFluor® 450.

REFERENCES

Willinger T and Flavell RA. 2012. Proc. Natl. Acad. Sci. 109:8670-8675. (flow cytometry) Meredith MM, Liu K, Darrasse-Jeze G, Kamphorst AO, Schreiber HA, Guernonprez P, Idoyaga J, Cheong C, Yao K-H, Niec RE, and Nussenzweig MC. 2012. J. Exp. Med. 209: 1153-1165. (immunofluorescence microscopy: acetone-fixed frozen tissue) Becker-Herman A, Meyer-Bahlburg A, Schwartz MA, Jackson SW, Hudkins KL, Liu C, Sather BD, Khim S, Liggitt D, Song W, Silverman GJ, Alpers CE and Rawlings DJ. 2011. J. Exp. Med. 208:2033-2042. (immunofluorescence microscopy – OCT embedded frozen tissue) Bertossi A, Aichinger M, Sansonetti P, Lech, Neff F, Pal M, Wunderlich FT, Anders H, Klein L, and Schmidt-Supprian M. 2011. J. Exp. Med. 208:1749-1756. (immunofluorescence microscopy) De Clercq S, Gembarska A, Denecker G, Maetens M, Naessens M, Haigh K, Haigh JJ, and Marine J-C. 2010. Mol. Cell. Biol. 30:5394-5405. (western blot) Nutt SL, Metcalf D, D'Amico A, Polli M, and Wu L. 2005. J. Exp. Med. 201:221-231. (immunomagnetic bead depletion) Cappione AJ, Pugh-Bernard AE, Anolik JH, and Sanz I. 2004. J. Immunol. 172: 4298-4307. (immunoprecipitation) Monteith CE, Chelack BJ, Davis WC, and Haines DM. 1996. Can. J. Vet. Res. 60(3): 193-198. (immunohistochemistry – feline tissue) Whiteland JL, Nicholls SM, Shimeld C, Easty DL, Williams NA, and Hill TJ. 1995. J. Histochem. Cytochem. 43:313-320. (immunohistochemistry – frozen tissue, paraffin embedded tissue) Domiati-Saad R, Ogle EW, and Justement LB. 1993. J. Immunol. 151: 5936-5947. (in vivo blocking)

NOTE: Please choose the appropriate format for each application. Citations are provided as a convenience to you; please consult Materials and Methods sections for additional details about the use of any product in these publications.

For Research Use Only.

Not for use in diagnostic or therapeutic procedures. Not for resale. Not for distribution without written consent. Tonbo Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Tonbo Biosciences, Tonbo Biosciences Logo and all other trademarks are the property of Tonbo Biotechnologies Corporation. © 2013 Tonbo Biosciences.