

TECHNICAL DATA SHEET

redFluor™ 710 Anti-Mouse CD4 (GK1.5)

Catalog Number: 80-0041

PRODUCT INFORMATION

Contents: redFluor™ 710 Anti-Mouse CD4 (GK1.5)

Isotype: Rat IgG2b, kappa

Concentration: 0.2 mg/mL


Clone: GK1.5

Reactivity: Mouse

Use By: 12 months from date of receipt

Storage Conditions: 2-8°C protected from light

Formulation: 10 mM NaH₂PO₄, 150 mM NaCl, 0.09% NaN₃, 0.1% gelatin, pH7.2

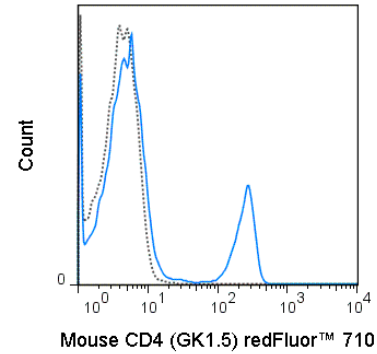

C57Bl/6 splenocytes were stained with 0.06 ug redFluor™ 710 Anti-Mouse CD4 (80-0041) (solid line) or 0.06 ug redFluor™ 710 Rat IgG2b (dashed line).

DESCRIPTION

The GK1.5 antibody reacts with mouse CD4, a 55 kDa protein which acts as a co-receptor for the T cell receptor (TCR) in its interaction with MHC Class II molecules on antigen-presenting cells. The extracellular domain of CD4 binds to the beta-2 domain of MHC Class II, while its cytoplasmic tail provides a binding site for the tyrosine kinase lck, facilitating the signaling cascade that initiates T cell activation. CD4 is typically expressed on thymocytes, certain mature T cell populations such as Th17 and T regulatory (Treg) cells, as well as on dendritic cells. The GK1.5 antibody is widely used as a phenotypic marker for CD4 expression. If used together, the GK1.5 antibody and an alternative antibody, Anti-Mouse CD4 clone RM4-5, will “compete” for binding, i.e. RM4-5 is able to block GK1.5 binding to cells. In contrast, the Anti-Mouse CD4 clone RM4-4 does not block binding of the GK1.5 antibody to cells (Arora S et al. 2006. Infect. Immun. 74: 4339-4348). The GK1.5 antibody is also reported to be cross-reactive with Syrian hamster CD4.

PREPARATION & STORAGE

This monoclonal antibody was purified from tissue culture supernatant via affinity chromatography. The purified antibody was conjugated under optimal conditions, with unreacted dye removed from the preparation. It is recommended to store the product undiluted at 4°C, and protected from prolonged exposure to light. Do not freeze.

APPLICATION NOTES

This antibody preparation has been quality-tested for flow cytometry using mouse spleen cells, or an appropriate cell type (where indicated). Please refer to the figure legend for the optimal concentration used to stain the tissue shown. We recommend titrating the antibody under your specific conditions to determine the optimal concentration of antibody needed in your experimental system.

redFluor™ 710 dye is excited by the red (633-647 nm) laser and has a peak emission of 710 nm. The recommended band pass filter for this dye is 710/50. redFluor™ 710 can be used as an alternative for Alexa Fluor® 700. Confirm that your cytometer is configured to detect this fluorochrome.

REFERENCES

Taniguchi RT, DeVoss JJ, Moon JJ, Sidney J, Sette A, Jenkins MK, and Anderson MS. 2012. Proc. Natl. Acad. Sci. 109: 7847-7852. (in vivo depletion)
 Wolkers MC, Gerlach, C, Arens R, Janssen EM, Fitzgerald P, Schumacher TN, Medema JP, Green DR, and Schoenberger SP. 2012. Blood. 119:798-804. (in vivo depletion)
 Lee L-F, Logronio K, Tu GH, Zhai W, Ni I, Mei L, Dilley J, Yu J, et al. 2012. Proc. Natl. Acad. Sci. 10:1073. (Flow cytometry).
 Thornton EE, Looney MR, Bose O, Sen D, Sheppard D, Locksley R, Huang X, and Krummel MF. 2012. J. Exp. Med. 10:1084. (Immunohistochemistry – OCT embedded frozen tissue)
 Stephen TL, Wilson BS, and Laufer TM. 2012. Proc. Natl. Acad. Sci. 109: 7415-7420. (Immunoprecipitation)
 Hammerbeck CD and Hooper JW. 2011. J. Virol. 85(19) 9929-9944. (Flow cytometry – Syrian hamster).
 Diamond MS, Kinder M, Matsushita H, Mashayekhi M, Dunn GP, Archambault JM, Lee H, Arthur CD, White JM, Kalinke U, Murphy KM, and Schreiber RD. 2011.

Tonbo Biosciences tests all antibodies by flow cytometry. Citations are provided as a resource for additional applications that have not been validated by Tonbo Biosciences. Please choose the appropriate format for each application and consult Materials and Methods sections for additional details about the use of any product in these publications.

For Research Use Only.
 Not for use in diagnostic or therapeutic procedures. Not for resale. Not for distribution without written consent. Tonbo Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Tonbo Biosciences, Tonbo Biosciences Logo and all other trademarks are the property of Tonbo Biotechnologies Corporation. © 2013 Tonbo Biosciences.