

TECHNICAL DATA SHEET

# redFluor™ 710 Anti-Mouse CD11c (N418)

Catalog Number: 80-0114

## PRODUCT INFORMATION

**Contents:** redFluor™ 710 Anti-Mouse CD11c (N418)

**Isotype:** Armenian Hamster IgG

**Concentration:** 0.2 mg/mL


**Clone:** N418

**Reactivity:** Mouse

**Use By:** 12 months from date of receipt

**Storage Conditions:** 2-8°C protected from light

**Formulation:** 10 mM NaH<sub>2</sub>PO<sub>4</sub>, 150 mM NaCl, 0.09% NaN<sub>3</sub>, 0.1% gelatin, pH7.2

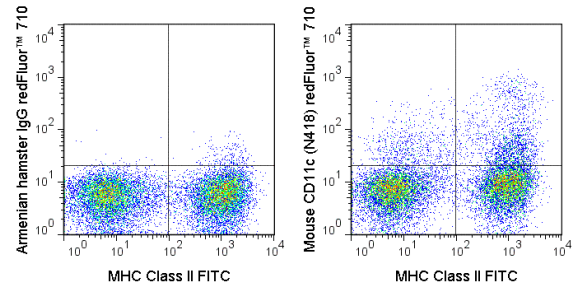

C57Bl/6 splenocytes were stained with FITC Anti-Mouse MHC Class II (35-5321) and 0.125 ug redFluor™ 710 Anti-Mouse CD11c (80-0114) (right panel) or 0.125 ug redFluor™ 710 Armenian Hamster IgG (left panel).

## DESCRIPTION

The N418 antibody reacts with mouse CD11c, also known as integrin alpha X. This 150 kDa cell surface glycoprotein is part of a family of integrin receptors that mediate adhesion between cells (cell-cell) and components of the extracellular matrix, e.g. fibrinogen (cell-matrix). In addition, integrins are active signaling receptors which recruit leukocytes to inflammatory sites and promote cell activation. Complete, functional integrin receptors consist of distinct combinations of integrin chains which are differentially expressed. Integrin alpha X (CD11c) assembles with Integrin beta-2 (CD18) into a receptor complex known as CR4 which can bind and induce signaling through ICAMs and VCAM-1 on endothelial cells and can also facilitate removal of iC3b bearing foreign cells. The N418 antibody is widely used as a marker for CD11c expression on dendritic cells (DC), often in parallel with markers for CD11b, for identification of developmental stages and mature subsets of this cell type. CD11c is prominently expressed on tissue macrophages, and is also detected on some types of activated T cells and intestinal intraepithelial lymphocytes (IEL).

## PREPARATION & STORAGE

This monoclonal antibody was purified from tissue culture supernatant via affinity chromatography. The purified antibody was conjugated under optimal conditions, with unreacted dye removed from the preparation. It is recommended to store the product undiluted at 4°C, and protected from prolonged exposure to light. Do not freeze.

## APPLICATION NOTES

This antibody preparation has been quality-tested for flow cytometry using mouse spleen cells, or an appropriate cell type (where indicated). The amount of antibody required for optimal staining of a cell sample should be determined empirically in your system.

redFluor™ 710 dye is excited by the red (633-647 nm) laser and has a peak emission of 710 nm. The recommended band pass filter for this dye is 710/50. redFluor™ 710 can be used as an alternative for Alexa Fluor® 700. Confirm that your cytometer is configured to detect this fluorochrome.

## REFERENCES

Guerrero JL, Ditsworth D, Catanzaro JM, Sabino G, Furie MB, Kew RR, Crawford HC, and Zong W-X. 2011. J. Immunol. 186: 3517-3526. (Immunohistochemistry – paraffin embedded tissue) Grewal JS, Pilgrim MJ, Grewal S, Kasman L, Werner P, Bruorton ME, London SD, and London L. 2011. FASEB J. 25:1680-1696. (Immunofluorescence microscopy – frozen tissue) Sadhu C, Ting HJ, Lipsky B, Hensley K, Garcia-Martinez LF, Simon SI, and Staunton DE. 2007. J. Leukoc. Biol. 81: 1395-1403. (in vitro blocking) Hagnerud, S, Manna PP, Cella M, Stenberg A, Frazier WA, Colonna M, and Oldenborg P-A. 2006. J. Immunol. 5772-5778. (Immunofluorescence microscopy – frozen tissue) Finkelman FD, Lees A, Birnbaum R, Gause WC, and Morris SC. 1996. J. Immunol. 157: 1406-1414. (in vivo activation) Huleatt JW and Lefrancois L. 1995. J. Immunol. 154: 5684-5693. (Immunoprecipitation) Metlay JP, Witmer-Pack MD, Agger R, Crowley MT, Lawless D, and Steinman RM. 1990. J. Exp. Med. 171: 1753. (Immunoprecipitation)

Tonbo Biosciences tests all antibodies by flow cytometry. Citations are provided as a resource for additional applications that have not been validated by Tonbo Biosciences. Please choose the appropriate format for each application and consult Materials and Methods sections for additional details about the use of any product in these publications.

For Research Use Only.

Not for use in diagnostic or therapeutic procedures. Not for resale. Not for distribution without written consent. Tonbo Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Tonbo Biosciences, Tonbo Biosciences Logo and all other trademarks are the property of Tonbo Biotechnologies Corporation. © 2013 Tonbo Biosciences.