

TECHNICAL DATA SHEET

violetFluor™ 500 Anti-Human CD3 (OKT3)

Catalog Number: 85-0037

PRODUCT INFORMATION

Contents: violetFluor™ 500 Anti-Human CD3 (OKT3)

Isotype: Mouse IgG2a, kappa

Concentration: 5 µL (0.5 µg)/test


Clone: OKT3

Reactivity: Human

Use By: 6 months from date of receipt

Storage Conditions: 2-8°C protected from light

Formulation: 10 mM NaH₂PO₄, 150 mM NaCl, 0.09% NaN₃,
0.1% gelatin, pH7.2

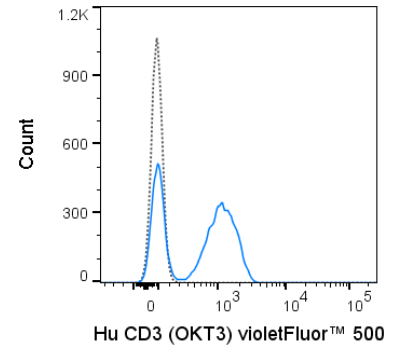

Human peripheral blood lymphocytes were stained with 5 µL (0.5 µg) violetFluor™ 500 Anti-Human CD3 (85-0037) (solid line) or 0.5 µg violetFluor™ 500 Mouse IgG2a isotype control (dashed line).

DESCRIPTION

The OKT3 antibody is specific for human CD3e, also known as CD3 epsilon, a 20 kDa subunit of the T cell receptor complex along with CD3 gamma and CD3 delta. These integral membrane protein chains assemble with additional chains of the T cell receptor (TCR), as well as CD3 zeta chain, to form the T cell receptor - CD3 complex. Together with co-receptors CD4 or CD8, the complex serves to recognize antigens bound to MHC molecules on antigen-presenting cells. These interactions promote T cell receptor signaling (T cell activation), inducing a number of cellular responses including proliferation, differentiation, production of cytokines or activation-induced cell death. CD3 is differentially expressed during thymocyte-to-T cell development and on all mature T cells. The OKT3 antibody is a widely used phenotypic marker for human T cells. In addition, as the CD3e subunit of the TCR complex contains intracellular signaling domains, binding of OKT3 antibody to CD3e can induce cell activation. A recent publication of the crystal structure of a CD3e-mitogenic antibody complex provides further insight as to the action of commonly used agonist antibodies, as well as specific epitope-binding data for the widely used human CD3 antibodies OKT3 and UCHT1 (Fernandes, R.A. et al. 2012. J. Biol. Chem. 287: 13324-13335). OKT3 has also been shown to be cross-reactive with Chimpanzee CD3 and

PREPARATION & STORAGE

This monoclonal antibody was purified from tissue culture supernatant via affinity chromatography. The purified antibody was conjugated under optimal conditions, with unreacted dye removed from the preparation. It is recommended to store the product undiluted at 4°C, and protected from prolonged exposure to light. Do not freeze.

APPLICATION NOTES

This antibody preparation has been pre-titrated and quality-tested for flow cytometry using an appropriate cell type. The antibody has been diluted for use at 5 µL per test, defined as the amount of antibody that will stain a cell sample in a final volume of approximately 100 µL. The number of cells within a sample should be determined empirically, but typically ranges between 1x10⁵ to 1x10⁸ cells.

violetFluor™ 500 dye is excited by the violet (405 nm) laser and has a peak emission of 500 nm. The recommended band pass filter for this dye is 525/20. violetFluor™ 500 can be used as an alternative for BD Horizon™ V500, eFluor® 506 or Pacific Orange®.

REFERENCES

- Haga-Friedman A, Horovitz-Fried M, and Cohen CJ. Jun 2012. J. Immunol. 188: 5538-5546. (in vitro activation)
- Bikker A, Moret FM, Kruize AA, Bijlsma JWJ, Lafeber FPJG, and van Roon JAG. Jun 2012. Rheumatology 51: 996-1005 (in vitro activation)
- Bagnara D, Kaufman MS, Calissano C, et al. 2011. Blood. 117: 5463-5472. (in vivo depletion)
- Nguyen V, Cao L, Lin JT, Hung N, Ritz A, Yu K, Jianu R, Ulin SP, Raphael BJ, Laidlaw DH, Brossay L, and Salomon AR. 2009. Mol. Cell. Proteomics. 8: 2418-2431. (in vitro activation)
- Bibollet-Ruche F, McKinney BA, Duverger A, Wagner FH, Ansari AA, and Kutsch O. 2008. J. Virol. 82(20): 10271-10278. (in vitro activation – Chimpanzee)
- Roura-Mir C, Catalfamo M, Cheng TY, Marqusee E, Besra GS, Jaraquemada D, and Moody DB. 2005. J. Immunol. 174: 3773-3780 (Immunohistochemistry – Acetone fixed frozen sections)

Tonbo Biosciences tests all antibodies by flow cytometry. Citations are provided as a resource for additional applications that have not been validated by Tonbo Biosciences. Please choose the appropriate format for each application and consult Materials and Methods sections for additional details about the use of any product in these publications.

For Research Use Only.

Not for use in diagnostic or therapeutic procedures. Not for resale. Not for distribution without written consent. Tonbo Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Tonbo Biosciences, Tonbo Biosciences Logo and all other trademarks are the property of Tonbo Biotechnologies Corporation. © 2013 Tonbo Biosciences.