

TECHNICAL DATA SHEET

violetFluor™ 500 Anti-Mouse CD4 (RM4-5)

Catalog Number: 85-0042

PRODUCT INFORMATION

Contents: violetFluor™ 500 Anti-Mouse CD4 (RM4-5)

Isotype: Rat IgG2a, kappa

Concentration: 0.2 mg/mL


Clone: RM4-5

Reactivity: Mouse

Use By: 12 months from date of receipt

Storage Conditions: 2-8°C protected from light

Formulation: 10 mM NaH₂PO₄, 150 mM NaCl, 0.09% NaN₃,
0.1% gelatin, pH7.2

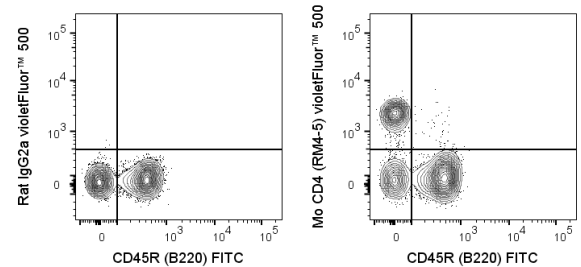

C57Bl/6 splenocytes were stained with FITC Anti-Mouse CD45R (B220) (35-0452) and 0.25 ug violetFluor™ 500 Anti-Mouse CD4 (85-0042) (right panel) or 0.25 ug violetFluor™ 500 Rat IgG2a isotype control (left panel).

DESCRIPTION

The RM4-5 antibody reacts with mouse CD4, a 55 kDa protein which acts as a co-receptor for the T cell receptor (TCR) in its interaction with MHC Class II molecules on antigen-presenting cells. The extracellular domain of CD4 binds to the beta2-domain of MHC Class II, while its cytoplasmic tail provides a binding site for the tyrosine kinase lck, facilitating the signaling cascade that initiates T cell activation. CD4 is typically expressed on thymocytes, certain mature T cell populations including Th17 and T regulatory (Treg) cells, as well as on dendritic cells. The RM4-5 antibody is widely used as a phenotypic marker for CD4 expression. If used together, the RM4-5 antibody and an alternative antibody, Anti-Mouse CD4 clone GK1.5, will compete for binding, i.e. RM4-5 antibody is able to block GK1.5 antibody binding to cells. In contrast, RM4-5 antibody does not block the binding of Anti-Mouse CD4 clone RM4-4 to cells.

PREPARATION & STORAGE

This monoclonal antibody was purified from tissue culture supernatant via affinity chromatography. The purified antibody was conjugated under optimal conditions, with unreacted dye removed from the preparation. It is recommended to store the product undiluted at 4°C, and protected from prolonged exposure to light. Do not freeze.

APPLICATION NOTES

This antibody preparation has been quality-tested for flow cytometry using mouse spleen cells, or an appropriate cell type (where indicated). Please refer to the figure legend for the optimal concentration used to stain the tissue shown. We recommend titrating the antibody under your specific conditions to determine the optimal concentration of antibody needed in your experimental system.

violetFluor™ 500 dye is excited by the violet (405 nm) laser and has a peak emission of 500 nm. The recommended band pass filter for this dye is 525/20. violetFluor™ 500 can be used as an alternative for BD Horizon™ V500 or Pacific Orange®.

REFERENCES

Willinger T and Flavell RA. 2012. Proc. Natl. Acad. Sci. 109:8670-8675. (Flow cytometry) Stephen TL, Wilson BS, and Laufer TM. 2012. Proc. Natl. Acad. Sci. 109: 7415-7420. (Immunofluorescence microscopy) Poitrasson-Riviere M, Biennu B, Le Campion A, Becourt C, Martin B, and Lucas B. 2008. J. Immunol. 180:7294-7304. (Immunohistochemistry – paraffin embedded tissue) Sorg H, Lorch B, Jaster R, Fitzner B, Ibrahim S, Holzhueter S, Nizze H, and Vollmar B. 2008. Am. J. Physio. Gastrointest. Liver Physiol. 295: G1274-1280. (Immunohistochemistry - paraffin embedded tissue) Menke J, Lucas JA, Zeller GC, Keir ME, Huang XR, Tsuboi N, Mayadas TN, Lan HY, Sharpe AH, and Kelley VR. 2007. J. Immunol. 179: 7466-7477. (Immunohistochemistry – frozen tissue) Irie J, Wu Y, Wicker LS, Rainbow D, Nalesnik MA, Hirsch R, Peterson LB, Leung PS, Cheng C, Mackay IR, Gershwin ME, and Ridgway WM. 2006. J Exp Med. 203(5):1209-19. (Immunohistochemistry – frozen tissue) Bosselut R, Zhang W, Ashe JM, Kopacz JL, Samelson LE, and Singer A. 1999. J. Exp. Med. 190: 1517-1526. (Immunoprecipitation) Shi Y, Kaliyaperumal A, Lu L, Southwood S, Sette A, Michaels MA, and Datta SK. 1998. J. Exp. Med. 187:367-378. (Blocking) Whiteland JL, Nicholls SM, Shimeld C, Easty DL, Williams NA, and Hill TJ. 1995. J. Histochem. Cytochem. 43:313-320. (Immunohistochemistry – frozen tissue, zinc-fixed paraffin embedded tissue)

Tonbo Biosciences tests all antibodies by flow cytometry. Citations are provided as a resource for additional applications that have not been validated by Tonbo Biosciences. Please choose the appropriate format for each application and consult Materials and Methods sections for additional details about the use of any product in these publications.

For Research Use Only.

Not for use in diagnostic or therapeutic procedures. Not for resale. Not for distribution without written consent. Tonbo Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Tonbo Biosciences, Tonbo Biosciences Logo and all other trademarks are the property of Tonbo Biotechnologies Corporation. © 2013 Tonbo Biosciences.